FOCUSLIGHT Never stop exploring

Laser Hair Removal

For cosmetic procedures where lasers are used to permanently remove unwanted hair, we offer optical components for beam shaping, as well as laser components and modules with integrated laser source, optics, and cooling.





Specifications	Value
Power Output	600 – 3000 W
Wavelength	755/808/1060 nm




Laser hair removal engines: Vsilk 2, Vsilk 2 Pro

Specifications	Vsilk 2	Vsilk 2 Pro
Power Output	Up to 1800 W	Up to 3000 W
Work Mode	106A, 30% Duty Cycle	185A, 10% Duty Cyc

Laser Body Sculpting

The non-invasive form of fat removal is performed by modules which consist of the laser source, and modules for beam shaping, cooling and skin contact detection.



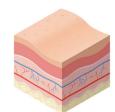

Skin and fat cells before and after laser body sculpting




Laser body sculpting module: Fairy

Fairy FR06 Features	Value
Power Output	50 W
Wavelength	1064 ± 20 nm
Beam Size	80 x 40 mm ²
Weight	385 g

Laser Skin Rejuvenation

Laser technology can significantly reduce skin imperfections like acne marks, scars, or wrinkles. Skin rejuvenation module with 1470nm laser is eye safe and therefore suitable for home-use application.



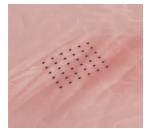

Laser skin rejuvenation module and the dot matrix it generated





FOCUSU 6HT

Fluorescence Angiography

Focuslight's fiber coupling modules and single emitter diode lasers can be applied in fluorescence endoscopy and other various laser medical treatment systems as core components with high reliability and excellent performance. A variety options of 808nm, 780nm, 650nm, and 450nm wavelength are provided.





Single-Emitter-Based Fiber Coupled Modules: FCMSE55

FCMSE55 Key Specs	Value
Power Output	15 - 25 W
Fiber Core Diameter	105/200 µm
Wavelength	808/780/650/450 nm
	·



Dust-Proof CW Single-Emitter Diode Laser: NV02


NV02 Key Specs	Value
Power Output	10 – 12 W
Fast Axis Collimation	1:1 Beam Ratio
Wavelength	808/780/650/450 nm

Disposable Endoscopy

Focuslight's high-quality wafer-level optics and wafer-level stacking technologies enable miniaturized micro-optical modules that can be easily integrated into disposable endoscopes. These precise elements improve the image quality and enhance the system performance.



Our Solution: Micro-Optical Modules Powered by Wafer-Level Optics and Stacking Technologies

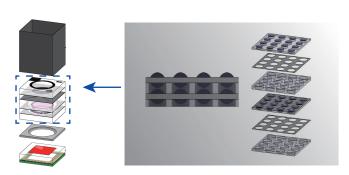

Key Features

- Tiny footprint of 1x1mm
- Full wafer-scale process for high volume production
- Micro-precision of wafer stacking
- Customization of optical performance available
- Leading thermal performance, ensuring simple thermal design


Health Monitoring

Focuslight optical components are used as an integral part of medical optical sensors, enabling real-time, non-invasive monitoring and device miniaturization in applications such as point of care, glucose monitoring and smart wearables.



Our Solutions for Health Monitoring:

- Shack-Hartmann Arrays
- Microlens Arrays
- Diffractive Optical Elements

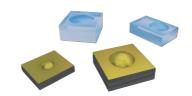






Confocal Microscopy


Focuslight produces high-precision lens discs and pinhole discs which are essential components of confocal microscopes. Our discs facilitate the generation of high-resolution, high-contrast, and sharply focused images for the 3D reconstruction of imaged samples.

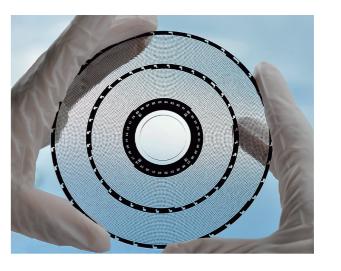

Pinhole Disc Features	Value
Disc diameter	≤ 160 mm
Material	Fused Silica
Pinhole / lens pattern	Custom
AR coating	UV, VIS, NIR
Chrome coating	Customized

Intraoral Scanner


Intraoral scanners typically employ structured light or laser projection systems to illuminate dental surfaces for 3D imaging. Diffractive optical elements (DOE) are used to shape or pattern the illumination, enabling detailed analysis of the dental structures. Microlens arrays (MLA) are integral to laser-based systems, where they help focus the laser beams for precise scanning.

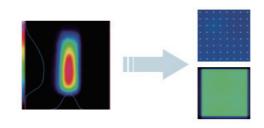

Spot generator and flat top beam





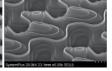

Diffractive optical elements

DOE Features	Value
Materials	FS and Si
Levels	2 (binary) to 16 levels
Wavelength Range	190nm to 5µm
Minimum Size	500nm to 1µm
Efficiency	Up to 96%

Our Solutions for Intraoral Scanners:

• Diffractive Optical Elements

• Microlens Array

MLA Features	Value
Angular Spectrum	Typical < 1 - 20 degrees
Materials	Fused Silica, Silicon
Area of Illumination	Linear, circular, rectangular, square
AR Coating	UV, VIS, NIR

Laser Surgery

Based on different absorption efficiencies of biological tissues, different wavelengths of laser are applied in various surgery applications. Focuslight offers various fiber coupled diode laser modules that create precise and concentrated laser output, which is essential in energy applications such as laser surgeries.



Urinary Surgery

Multi-Bar-Based Fiber

Power Output 100 - 400 W (QCW)

808/980 nm

Coupled Modules





Side-Pumped Modules

Power Output



Single-Emitter-Based

Fiber Coupled Modules

Dentistry





Single-Bar-Based Fiber Coupled Modules

Eye Surgery

Specifications	Value
Power Output	25 - 35 W
Fiber Core Diameter	200 µm
Wavelength	808/976/ 1064 nm

	Specifications	Value
	Power Output	10 – 12 W
	Fast Axis Collimation	1:1 Beam F Available
	Wavelength	808/915/9 976 nm

Open Package Single

Emitter Diode Laser

Specifications	Value
Power Output	15 - 25 W
Fiber Core Diameter	400 µm
Wavelength	808 nm